
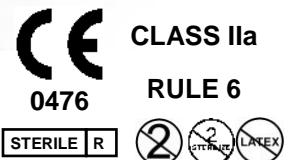


Ball tipped suction cannula for NEUROSURGERY

Complete with ergonomic aspiration control handle

- Aspiration control hole has a diameter of 1,5 mm., allowing the surgeon a slower aspiration response. Whenever the hole is opened the aspiration on the tip will decrease progressively over time.
- Particularly adapt for Neurosurgery operations.
- The ball shaped distal tips reduce tissue trauma.
- Made of stainless medical grade AISI 304 steel.
- Angled, with uniform internal diameter over the entire cannula length.
- Handle with aspiration control made of copolymer styrene-butadiene (ABS).
- Color coded for immediate recognition.
- Latex free.
- Sterilized by gamma rays.
- Packaged in double sterile envelopes.

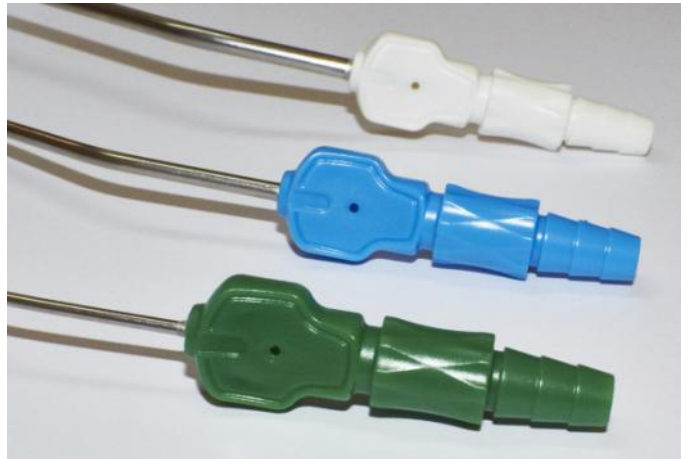

Code	Internal diameter mm.	External diameter mm.	CH	Cannula length, mm.	Handle colour	Pcs per BOX
CSSR20Y130*	1,45	2,0	6	130	GREEN	240
CSSR20Y160*	1,45	2,0	6	160	GREEN	240
CSSR30Y130*	2,54	3,0	9	130	BLUE	240
CSSR30Y160*	2,54	3,0	9	160	BLUE	240
CSSR40Y130*	3,0	4,0	12	130	WHITE	240
CSSR40Y160*	3,0	4,0	12	160	WHITE	240

* Product not in stock, available by specific request

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
 Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
 www.simitalia.it www.medeurope.eu
 E-Mail: info@medeurope.eu

Distributed by:



Ball tipped suction cannula for NEUROSURGERY

Complete with ergonomic drop shaped aspiration control handle



CLASS IIa

0476

RULE 6

STERILE R



- The aspiration control hole has an ergonomic drop shaped form which allows the surgeon high precision control of the aspiration pressure. When the hole is opened, the aspiration at the tip decreases progressively over time.
- Particularly adapt for Neurosurgery operations
- The ball shaped distal tips reduce tissue trauma
- Cannulae are made of stainless medical grade AISI 304 steel
- Uniform internal diameter over the entire cannula length
- Angled
- Handle with aspiration control made of copolymer styrene-butadiene (ABS)
- Color coded for immediate recognition
- Does not contain latex
- Sterilized by Gamma Rays
- Packaged in double sterile envelopes

Code	Internal Diameter mms.	External Diameter mms.	CH	Total length Of Cannula mms.	Handle Colour	Pcs. per BOX
CSSR20B130*	1,45	2,0	6	130	GREEN	240
CSSR30B130*	2,54	3,0	9	130	BLUE	240
CSSR40B130*	3,0	4,0	12	130	WHITE	240
CSSR20B160*	1,45	2,0	6	160	GREEN	240
CSSR30B160*	2,54	3,0	9	160	BLUE	240
CSSR40B160*	3,0	4,0	12	160	WHITE	240

*Product not in stock, available on request

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
 Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
 www.simalia.it www.medeurope.eu
 E-Mail: info@medeurope.eu

Distributed by:

