
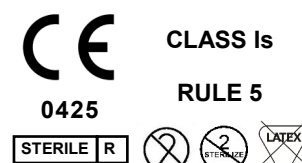


Med Stomach Tube

Gastric calibration tube

Revision date: 31.08.2026

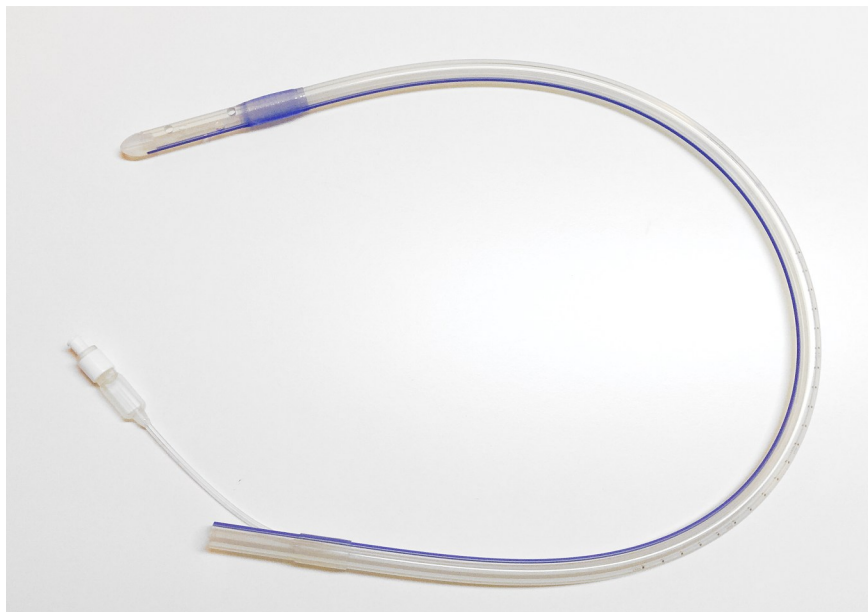

The device is indicated for use in surgical procedures unique to the stomach. Introduced by mouth of the patient it allows, through its asymmetrical distal balloon, the movement of tissues, facilitating the surgeon involved in the intervention in laparoscopic or open surgery. The holes on the distal portion of atraumatic tip favor the complete suction.

Characteristics

- Material: medical grade silicone
- With radiopaque line
- With or without asymmetrical distal balloon
- Length 800 mm. or 1000mm
- Sterilised by gamma rays
- Latex and Phthalates free
- Packaging: outer sterile pouch and inner protective bag.



Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850
www.simitalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:

arfi

CODE / REF	CH	Shore	Length	Model
MST-36-AC-01-080*	36	70	80 cm	8 distal side holes Ø 5 mm
MST-38-AC-01-080*	38	70	80 cm	8 distal side holes Ø 5 mm
MST-34-AF-01-080*	34	80	80 cm	8 distal side holes Ø 5 mm
MST-36-AF-01-080*	36	80	80 cm	8 distal side holes Ø 5 mm
MST-38-AF-01-080*	38	80	80 cm	8 distal side holes Ø 5 mm
MST-40-AF-01-080*	40	80	80 cm	8 distal side holes Ø 5 mm
MST-42-AF-01-080*	42	80	80 cm	8 distal side holes Ø 5 mm
MST-44-AF-01-080*	44	80	80 cm	8 distal side holes Ø 5 mm
MST-36-BC-01-080*	36	70	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-38-BC-01-080*	38	70	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-34-BL-01-080*	34	80	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-36-BL-01-080*	36	80	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-38-BL-01-080*	38	80	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-40-BL-01-080	40	80	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-42-BL-01-080*	42	80	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-44-BL-01-080*	44	80	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-36-BN-01-080*	36	57	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-38-BN-01-080*	38	57	80 cm	with asymmetrical distal balloon and 3 distal side holes Ø 5 mm
MST-34-LL-01-080*	34	57	80 cm	5 distal side holes Ø 5 mm
MST-36-LL-01-080*	36	57	80 cm	5 distal side holes Ø 5 mm
MST-38-LL-01-080*	38	57	80 cm	5 distal side holes Ø 5 mm
MST-40-LL-01-080*	40	57	80 cm	5 distal side holes Ø 5 mm
MST-42-LL-01-080*	42	57	80 cm	5 distal side holes Ø 5 mm
MST-44-LL-01-080*	44	57	80 cm	5 distal side holes Ø 5 mm
MST-34-LL-02-100*	34	57	100 cm	with black tip and 5 distal side holes Ø 5 mm
MST-36-LL-02-100*	36	57	100 cm	with black tip and 5 distal side holes Ø 5 mm
MST-38-LL-02-100*	38	57	100 cm	with black tip and 5 distal side holes Ø 5 mm
MST-40-LL-02-100*	40	57	100 cm	with black tip and 5 distal side holes Ø 5 mm
MST-42-LL-02-100*	42	57	100 cm	with black tip and 5 distal side holes Ø 5 mm
MST-44-LL-02-100*	44	57	100 cm	with black tip and 5 distal side holes Ø 5 mm

* Product not in stock, available by specific request

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
 Tel. +39 051.817902 +39 051.817850
 www.simitalia.it www.medeurope.eu
 E-Mail: info@medeurope.eu

Distributed by:

