



NEW Safety Count System


maxi box for special surgery

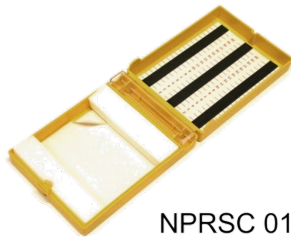


NPRSC 013

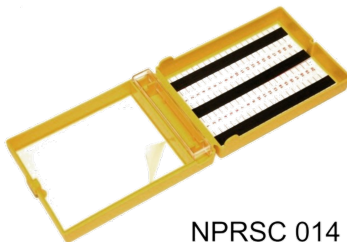


NPRSC 014

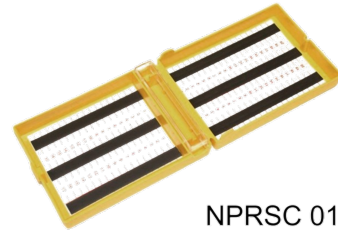


NPRSC 015

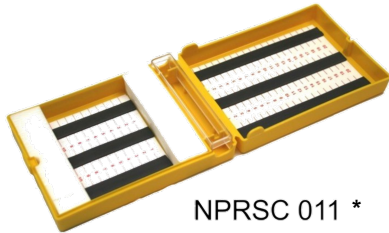


NPRSC 011 *

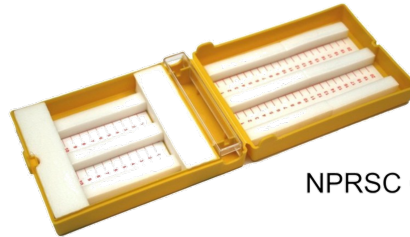


NPRSC 017 *

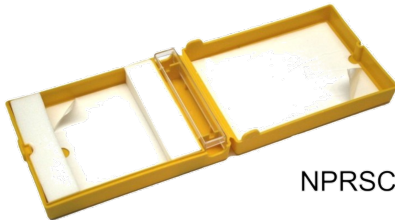


NPRSC 012 *

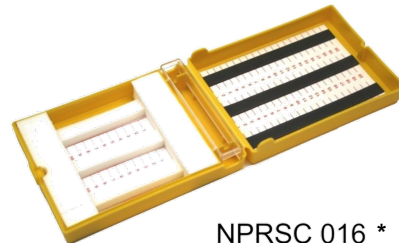

NPRSC 016 *

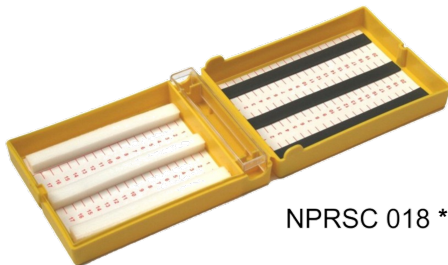

NPRSC 018 *

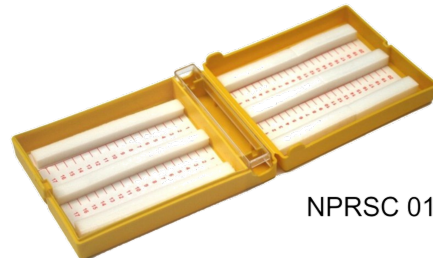

NPRSC 019 *



CLASS Is

RULE 1

STERILE R



* Product not in stock, available by specific request

Indications

The New Safety Count System Maxi Box series permits the highest standard of quality elimination of scalpel blades and suture needles used in operating rooms under the current European standards of the Medical Devices.

Characteristics

- The New Safety Count System Maxi series permits quality counting of scalpel blades and needles used in the operating room, and offers the highest safety standards to avoid the risk of contamination to dispose of sharps of all types.
- They are safe, simple, easy to use and economic
- Various models are available to offer convenience of choice to the user
- Single sterile packaging in medical grade paper and polypropylene (sterilised by gamma rays)
- Dimensions open: 26,5 cms. x 9,5 cms. x 2,5 cms.
- Dimensions closed: 13,2 cms. x 9,5 cms. x 3,5 cms.
- Box of 60 pcs.

The blade detachment device found on either side of the device makes use simple and safe for both left handed and right handed users

Numbered background for precise counting

Transparent polystyrene cover

Foam water based adhesive bottom that retains contaminated materials

Magnetic strips or foam water based adhesive to retain used suture needles

After inserting the scalpel into the blade detachment system, push down and extract the handle

Safety click closure that resists accidental opening but which may be re-opened using the scalpel handle

Extended internal border for increased safety

Resistent container made of yellow, puncture safe polypropylene compound

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850
www.simitalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:

