
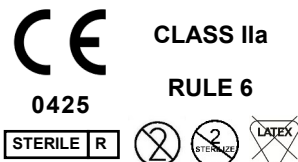


Ball tipped suction cannula

Complete with ergonomic aspiration control handle

- Aspiration control hole has a diameter of 1,5 mm., allowing the surgeon a slower aspiration response. Whenever the hole is opened the aspiration on the tip will decrease progressively over time.
- The ball shaped distal tips reduce tissue trauma.
- Made of stainless medical grade AISI 304 steel.
- Angled, with uniform internal diameter over the entire cannula length.
- Suction control handle made of SBC plastic material
- Color coded for immediate recognition.
- Latex free.
- Sterilized by gamma rays.
- Packaged in double sterile envelopes.

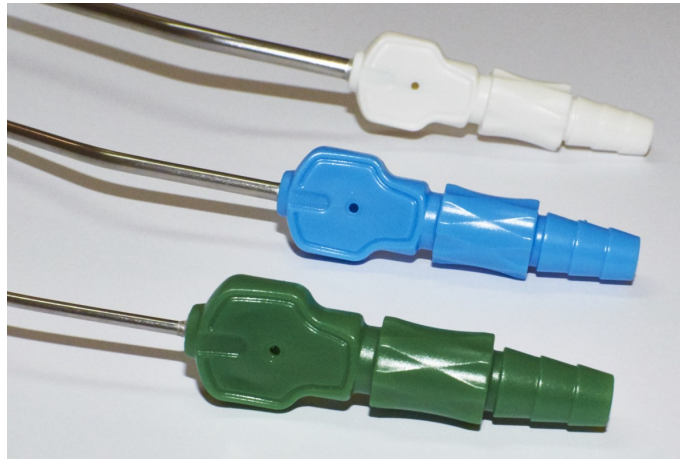

| Code | Internal diameter mm. | External diameter mm. | CH | Cannula length, mm. | Handle colour | Pcs per BOX |
|-------------|-----------------------|-----------------------|----|---------------------|---------------|-------------|
| CSSR20Y130* | 1,45 | 2,0 | 6 | 130 | GREEN | 240 |
| CSSR20Y160* | 1,45 | 2,0 | 6 | 160 | GREEN | 240 |
| CSSR30Y130* | 2,54 | 3,0 | 9 | 130 | BLUE | 240 |
| CSSR30Y160* | 2,54 | 3,0 | 9 | 160 | BLUE | 240 |
| CSSR40Y130* | 3,0 | 4,0 | 12 | 130 | WHITE | 240 |
| CSSR40Y160* | 3,0 | 4,0 | 12 | 160 | WHITE | 240 |

* Product not in stock, available by specific request

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
www.simitalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:



Ball tipped suction cannula

Complete with ergonomic drop shaped aspiration control handle



- The aspiration control hole has an ergonomic drop shaped form which allows the surgeon high precision control of the aspiration pressure. When the hole is opened, the aspiration at the tip decreases progressively over time.
- The ball shaped distal tips reduce tissue trauma
- Cannulae are made of stainless medical grade AISI 304 steel
- Uniform internal diameter over the entire cannula length
- Angled
- Suction control handle made of SBC plastic material
- Color coded for immediate recognition
- Does not contain latex
- Sterilized by Gamma Rays
- Packaged in double sterile envelopes



CLASS IIa

0425

RULE 6

STERILE R



| Code | Internal Diameter mms. | External Diameter mms. | CH | Total length Of Cannula mms. | Handle Colour | Pcs. per BOX |
|-------------|------------------------|------------------------|----|------------------------------|---------------|--------------|
| CSSR20B130* | 1,45 | 2,0 | 6 | 130 | GREEN | 240 |
| CSSR30B130* | 2,54 | 3,0 | 9 | 130 | BLUE | 240 |
| CSSR40B130* | 3,0 | 4,0 | 12 | 130 | WHITE | 240 |
| CSSR20B160* | 1,45 | 2,0 | 6 | 160 | GREEN | 240 |
| CSSR30B160* | 2,54 | 3,0 | 9 | 160 | BLUE | 240 |
| CSSR40B160* | 3,0 | 4,0 | 12 | 160 | WHITE | 240 |

*Product not in stock, available on request

Manufactured by:




Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
www.simalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:



Ball Tipped Aspiration Cannula for ENT

Complete with suction control handle

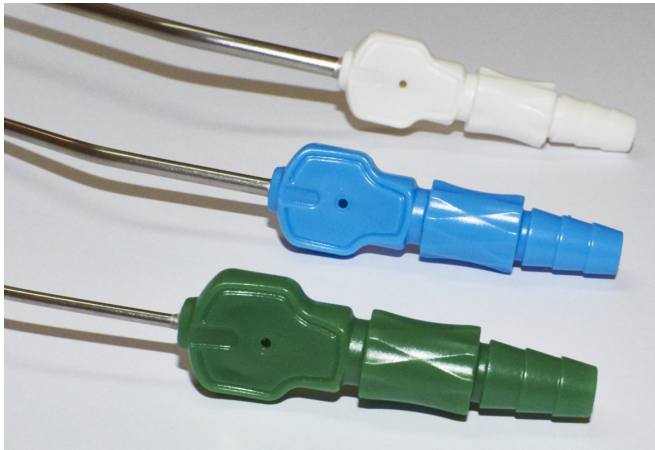

- Aspiration control hole has a diameter of 1,5 mm., allowing the surgeon a slower aspiration response. Whenever the hole is opened the aspiration on the tip will slowly decrease. Even with a very high aspiration pressure the cannula makes only a low amount of noise.
- Particularly adapt for ENT applications
- The ball shaped distal tips reduce tissue trauma
- Cannula are made of stainless medical grade AISI 304 steel
- Uniform internal diameter over the entire cannula length
- Angled
- Suction control handle made of SBC plastic material
- Color coded for immediate recognition
- Does not contain latex
- Sterilized by Gamma Rays
- Packaged in single sterile envelopes



| Code | Internal Diameter mms. | External Diameter mms. | CH | Total length Of Cannula mms. | Handle Colour | Pcs. per BOX | Pcs. per EXPORT BOX |
|-------------|------------------------|------------------------|----|------------------------------|---------------|--------------|---------------------|
| CSSR200130* | 1,45 | 2,0 | 6 | 130 | GREEN | 60 | 240 |
| CSSR300130* | 2,54 | 3,0 | 9 | 130 | BLUE | 60 | 240 |
| CSSR400130* | 3,0 | 4,0 | 12 | 130 | WHITE | 60 | 240 |
| CSSR200160* | 1,45 | 2,0 | 6 | 160 | GREEN | 60 | 240 |
| CSSR300160* | 2,54 | 3,0 | 9 | 160 | BLUE | 60 | 240 |
| CSSR400160* | 3,0 | 4,0 | 12 | 160 | WHITE | 60 | 240 |

*Product not in stock, available on request

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
www.simitalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:

