


NEW Safety Count System

Revision date: 18.04.2023

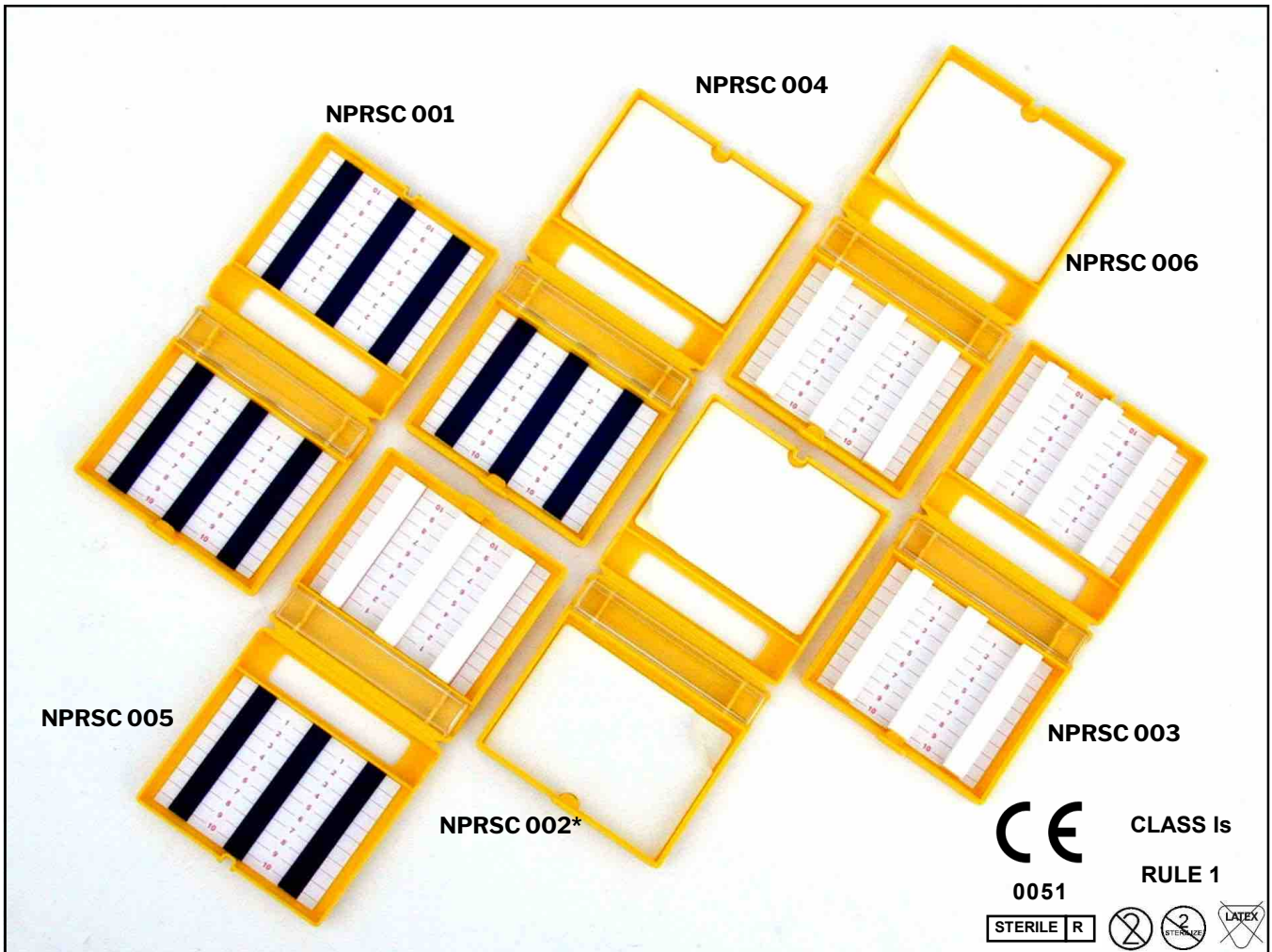

The SAFETY COUNT SYSTEM sharps counting and disposal system guarantees the safe elimination of scalpel blades and suture needles used in the operating room.

Characteristics

Used sharps of all kinds can be counted and disposed of without the risk of contamination. The Safety Count System is safe, simple and economic.

- different models allow maximum user choice;
- the exclusive bilateral blade detachment system is efficient on a large variety of scalpel blades
- incorporates a resistant safety closure which protects against accidental opening (which can however be opened if necessary using the scalpel handle) ;
- strong, yellow, plastic "anti-puncture" container;
- sterile, blister packaged;
- conveniently packaged 18 pcs./dispenser box;
- sterilized by gamma rays.







Codes

- NPRSC 001: magnet + magnet
- NPRSC 002*: adhesive + adhesive
- NPRSC 003: foam + foam
- NPRSC 004: adhesive + magnet
- NPRSC 005: foam + magnet
- NPRSC 006: foam + adhesive

Dispenser containing 18 Safety Count System units.
108 pcs. per box

* Product not in stock, available by specific request

PATENT PENDING

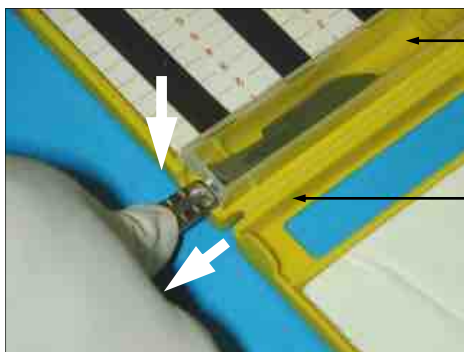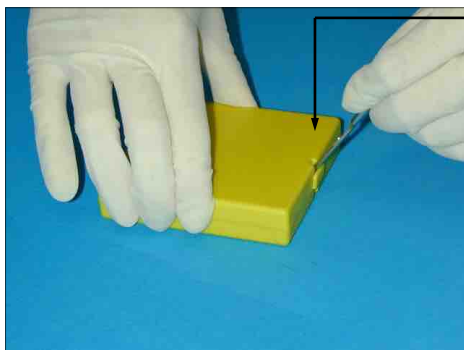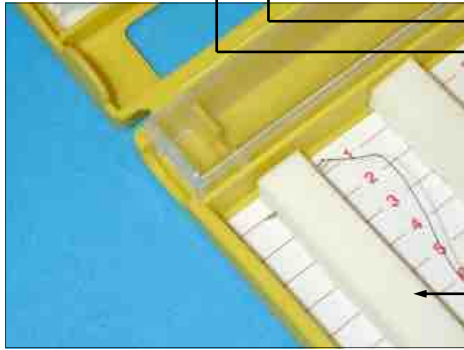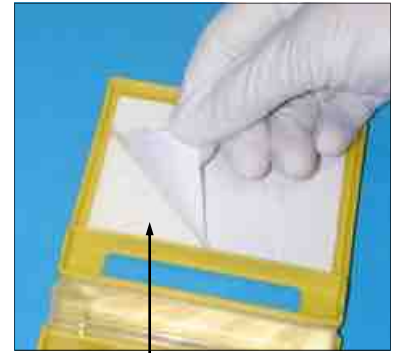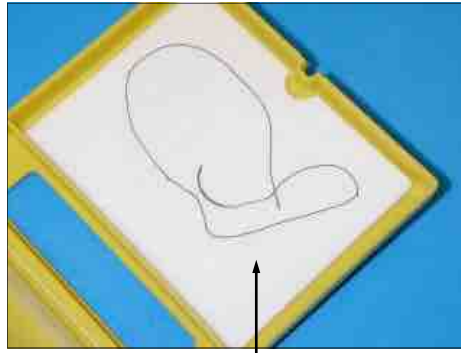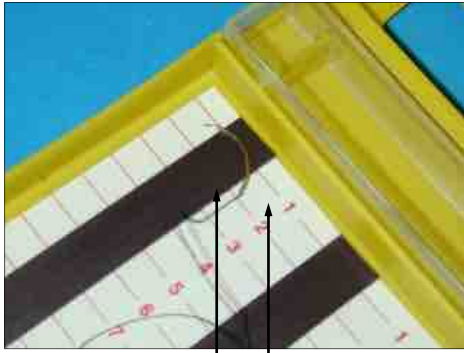

- Adhesive base which withholds contaminated materials
- Numbered base for precise counting
- Magnetic or foam strips that withhold used suture needles
- Strong, punctureproof yellow plastic container
- Safe click-locking closure which resists accidental opening, but which may be opened intentionally using the scalpel handle
- The double ended blade detachment device, provides easy and safe use for both right and left handed users
- Transparent cover
- After inserting the scalpel into the detachment device, push and pull back on the handle
- Self-fitting internal edges for increased closure safety
- Practical 18 piece dispenser

CE
0051

STERILE R



CLASS Is
RULE 1

PATENT PENDING



Manufactured by:




Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
www.simitalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:




NEW Safety Count System


maxi box for special surgery

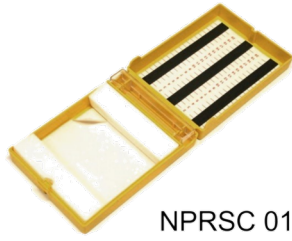


NPRSC 013

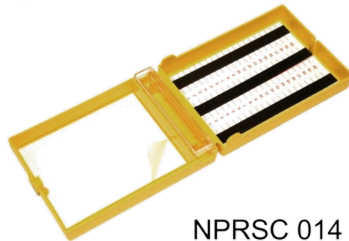


NPRSC 014

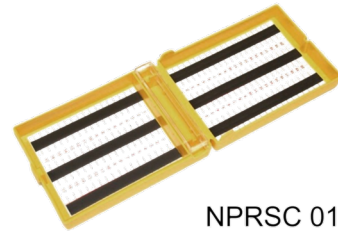


NPRSC 015

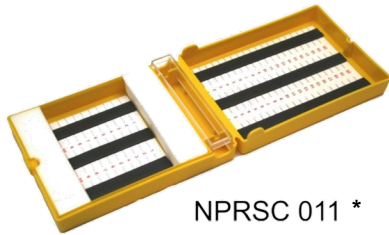


NPRSC 011 *

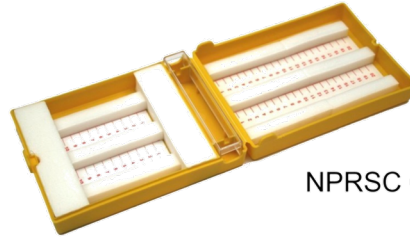


NPRSC 017 *

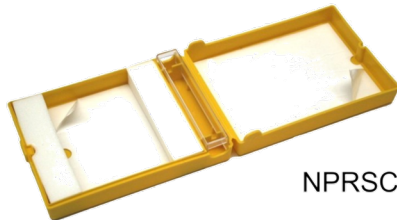


NPRSC 012 *

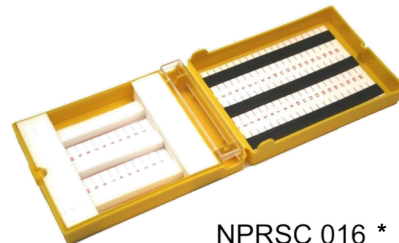

NPRSC 016 *

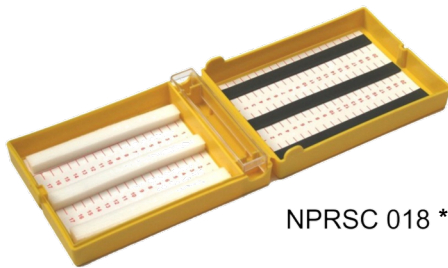


NPRSC 018 *

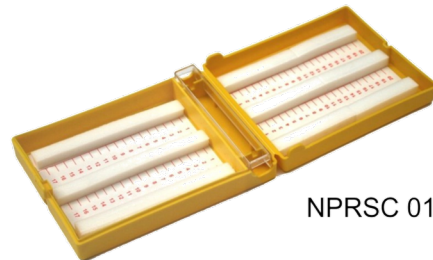

NPRSC 019 *

* Product not in stock, available by specific request

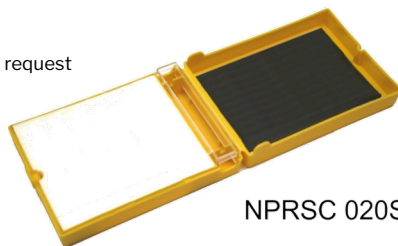


NPRSC 020S *



CLASS Is

RULE 1

STERILE R

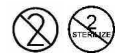

Indications

The New Safety Count System Maxi Box series permits the highest standard of quality elimination of scalpel blades and suture needles used in operating rooms under the current European standards of the Medical Devices.

Characteristics


- The New Safety Count System Maxi series permits quality counting of scalpel blades and needles used in the operating room, and offers the highest safety standards to avoid the risk of contamination to dispose of sharps of all types.
- They are safe, simple, easy to use and economic
- Various models are available to offer convenience of choice to the user
- Single sterile packaging in medical grade paper and polypropylene (sterilised by gamma rays)
- Dimensions open: 26,5 cms. x 9,5 cms. x 2,5 cms.
- Dimensions closed: 13,2 cms. x 9,5 cms. x 3,5 cms.
- Box of 60 pcs.

The blade detachment device found on either side of the device makes use simple and safe for both left handed and right handed users

Numbered background for precise counting

Transparent polystyrene cover


Foam water based adhesive bottom that retains contaminated materials

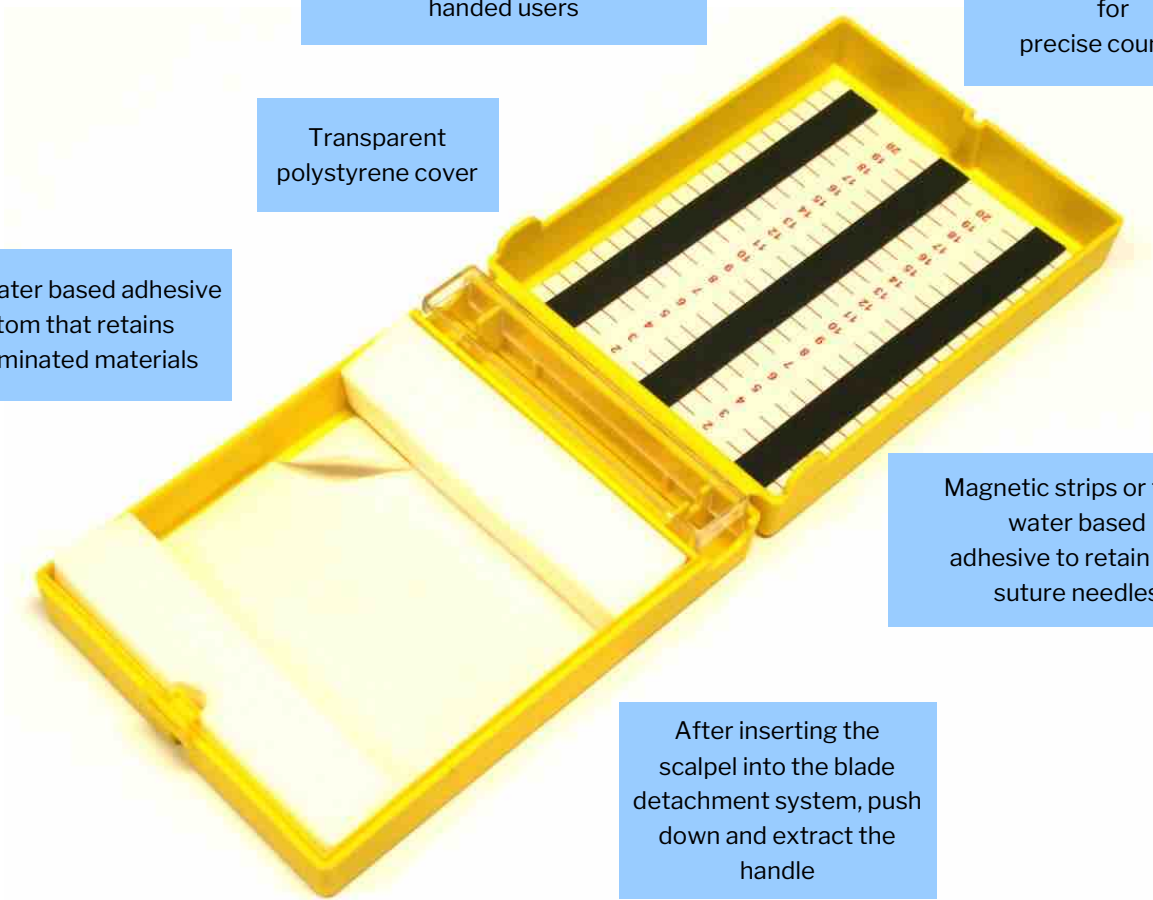

Magnetic strips or foam water based adhesive to retain used suture needles

After inserting the scalpel into the blade detachment system, push down and extract the handle

Safety click closure that resists accidental opening but which may be re-opened using the scalpel handle

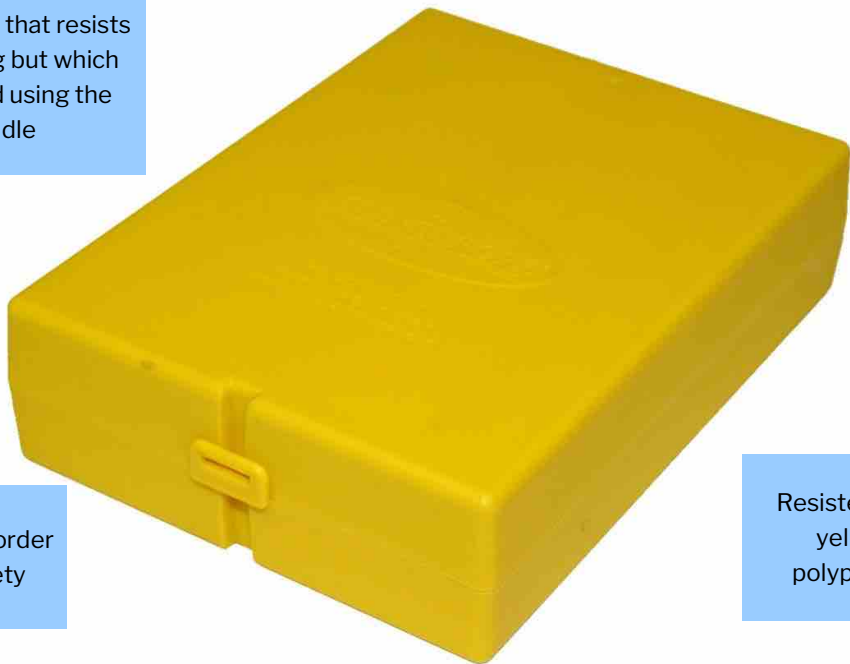

Resistent container made of yellow, puncture safe polypropylene compound

Extended internal border for increased safety

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850 Fax +39 051.817845
www.simitalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:

