


GASTRIC BAND

AllBand

Ring For Bariatric Surgery

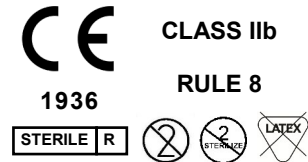

Implantable silicone ring with adjustable circumference used in bariatric surgery procedures, such as banded sleeve gastrectomy, gastric bypass, and banded mini-gastric bypass. The device features a graduated scale and a black pad-printed side of the tip (distal portion) to assist the surgeon in inserting and positioning the ring around the stomach portion. ALLBAND can be implanted using both laparoscopic and laparotomic techniques. The device is indicated for both primary and revision surgery. The device must be secured with non-absorbable stitches to prevent slippage.

Characteristics


- Section measures: 2.5x4 mm.
- Length: 205 mm.
- Circumference corresponding to internal diameter: see the table.
- Material: medical grade implantable silicone, white color, x-ray, shore: 70A.
- Sterilised by gamma rays.
- Double packaging (internal blister and external pouch).
- Latex free.

CIRCUMF. mm	INTERNAL DIAM. mm
48,7	15
63,7	20
78,7	25
93,7	30
108,7	35
123,7	40

Code

REF **BB 13-45 P**

Packaging: 1 piece/box

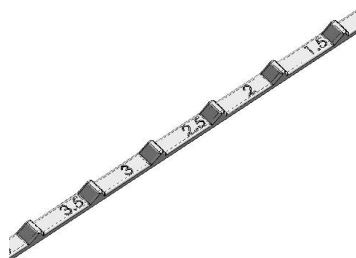

Grooved section for device closure

Manufactured by:



Via Gramsci, 9 - 40018 San Pietro in Casale (BO) - Italy
Tel. +39 051.817902 +39 051.817850
www.simitalia.it www.medeurope.eu
E-Mail: info@medeurope.eu

Distributed by:

